


CE

Aegistab OP850 Operating Table

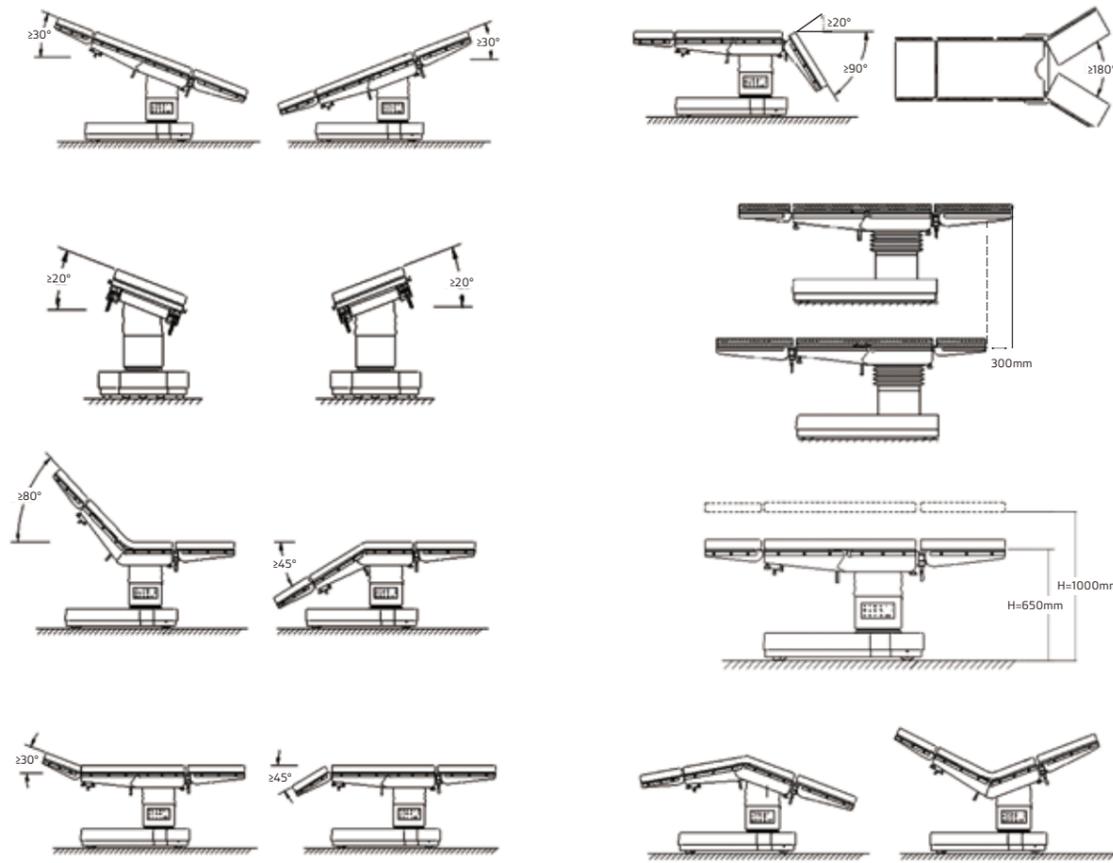

Technical Specifications

Model	OP850
Table Length (mm)	2080
Table Width (w/o side rails, mm)	520
Table Height (w/o pads, mm)	650 - 1000
Trendelenburg/Reverse	30°/30°
Lateral Tilt (left/right)	20°/20°
Back Plate Up/Down	80°/45°
Head Plate Up/Down (detachable)	30°/45°
Leg Plate Up/Down (detachable)	20°/90°
Leg Plate Split (detachable)	180°
Flex/Reflex	225° / 110°
Kidney-bridge Up (mm)	128
Longitudinal Shift (mm)	300 (Optional)
Weight Capacity (kg)	350
Exchangeable Head and Leg Plates	Yes
One Key Reset	Yes





Corporate Headquarter
HEYER Medical AG
 Carl-Heyer-Str. 1/3, 56130 Bad Emms, Germany
 Tel: +49 2603 791 0
 Fax: +49 2603 791 200
 info@heyemed.com
 www.heyemedical.de

Manufacturer
HEYER Medical Co., Ltd.
 No.10, Chaobai Street, Yanjiao Development Zone
 065201, Sanhe City, Hebei, China
 Tel: +86 10 5841 1198
 Fax: +86 10 6371 8989
 marketing@heyemed.com



An excellent electro-hydraulic operating table that can support up to 350 kg of weight and has multiple surgical configurations to accommodate various surgical needs.

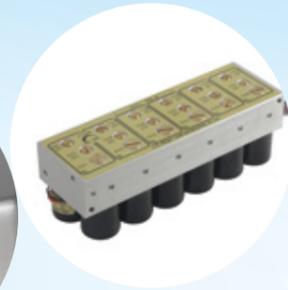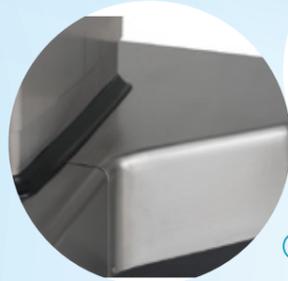

① Premium Quality Material

- All exposed parts are made of HI-Q 304 stainless steel
- High corrosion resistance, easy to clean and disinfection



② Built-in Elevator Kidney-bridge

- Built-in X-ray radiolucent kidney-bridge supports kidney and gall bladder surgeries efficiently



③ Perfect X-Ray Imaging

- Head and leg plates with excellent radiolucent guarantee high quality X-ray imaging
- Exchangeable head & leg sections
- Maximum 300 mm longitudinal shifting function (optional)



④ Surgical Configuration

- One-key reset function
- Flex and reflex function
- Flexible patient positionings with optional accessories are suitable for a wide variety of surgeries



⑤ Ergonomic Adjustment

- Convenient leg plate and head plate positionings ensure one-hand control operation
- Imported gas springs guarantee smooth movement and accurate positioning
- Wide range of table height adjustment and T-shape base provide comfortable working environment



⑥ Safety

- Auto-lock system: table base will be locked automatically once the position was adjusted by controller, preventing possible misoperation
- Two control systems: remote hand controller and column control panel work independently, if one is broken, another still work as normal
- Intelligent safety system: Control panel will power off automatically when there is no operation over 1 minute, avoiding any possible misoperation



⑦ Pads

- High density memory foam with soft texture can effectively reduce the occurrence of bedsores